



DUOMED



G-Lock

DUODENOSCOPE ED-580XT

FUJIFILM
Value from Innovation




GASTROENTEROLOGY

VIDEO DUODENOSCOPE **ED-580XT**

ERCP (Endoscopic retrograde cholangiopancreatography) uses a combination of endoscopy and fluoroscopy to diagnose and treat the biliary or pancreatic ductal system and requires a high level of expertise from the endoscopist.

The Fujifilm video duodenoscope is designed to support quality and efficiency and utilises a new G-Lock elevator. A unique distal attachment cap allows easy, fast and efficient cleaning. The ED-580XT is made to support you during your most complex procedures.

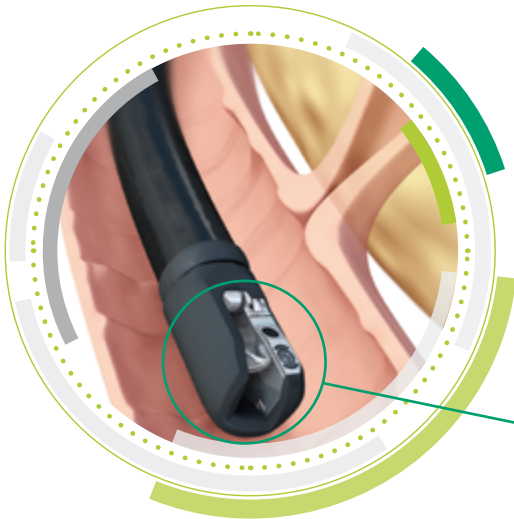

G-Lock

Improved treatment capability

**ADVANCED FORCE
TRANSMISSION &
MANOEUVRABILITY**

**GREATER
ELEVATION
FORCE**

**EASY AND
EFFECTIVE
DISTAL END
CLEANING**





IMPROVED WIRE LOCKING **G-LOCK** & GREATER **ELEVATION FORCE**



IMPROVED TREATMENT CAPABILITY

Incorporated into the distal tip of the ED-580XT, the G-Lock contains the forceps elevator and the contact section, enabling the guidewire to be simply and securely fixed into position by using the forceps elevator. In addition, the specially designed round shaped forceps elevator reduces the risk of guidewire damage. The inner tube of the instrument channels uses an optimised material to enable a device to be inserted smoothly, supporting rapid device exchange. Designed to work in harmony with the endoscopist, the new G-Lock and low friction instrument channel support efficiency and ease of use during ERCP procedures.

INCREASED ANGLE OF ENDOTHERAPY DEVICES

The forceps elevator design increases the maximum angle that endotherapy devices can be introduced. Movement of the elevator lever at the control body of the duodenoscope translates accurately to the distal tip, allowing for very precise positioning of instruments.




Case image of improved elevation force

ADVANCED MANOEUVRABILITY

ADVANCED FORCE TRANSMISSION




The flexible portion is designed to transmit the pushing, pulling and rotating movements from the hand to the distal end of the endoscope, providing enhanced manoeuvrability inside the digestive tract.

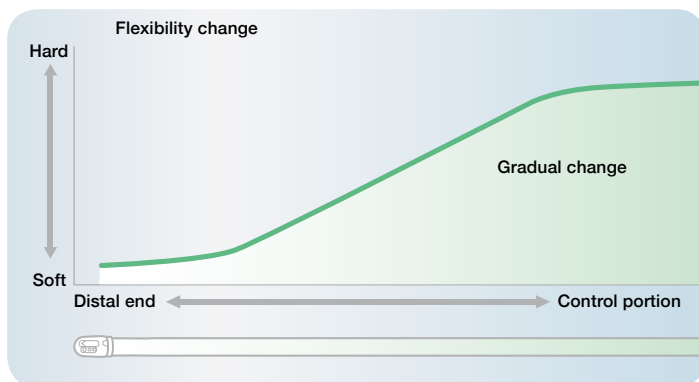

INTUITIVE OPERATION



The G7 control portion has a rounded surface design to fit in the hand, and its layout makes intuitive operation possible.

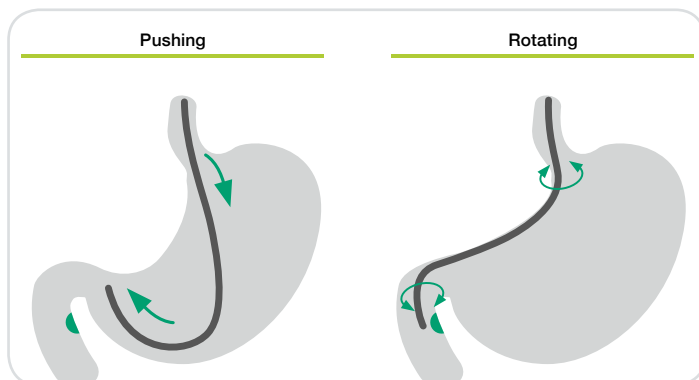



EASY AND EFFECTIVE DISTAL END CLEANING



EASIER BRUSHING ACCESS – EASIER CLEANING

The single-use distal end cap permits easier brushing access to the distal end of the endoscope. In addition, the elevator mechanism is sealed to allow easier cleaning.

ATTACHMENT AND DETACHMENT OF DISTAL END CAP

ATTACHMENT



Cover the distal end cap until the installation groove of the distal end is hidden completely.



Make sure that the distal end cap is attached tightly by pulling it straight.

DETACHMENT



Lift up a portion of distal end cap near the instrument channel outlet.




While turning the portion counterclockwise, pull it straight out.

SPECIFICATIONS

ED-580XT DUODENOSCOPE



Viewing direction	95° (retro-viewing 5°)
Field of view	100°
Observation range	4–60 mm
Distal end diameter	13.1 mm
Bending capability	Up 120° / Down 90° Left 90° / Right 110°
Working length	1,250 mm
Total length	1,550 mm
Insertion tube diameter	11.3 mm
Minimum diameter of instrument channel	4.2 mm

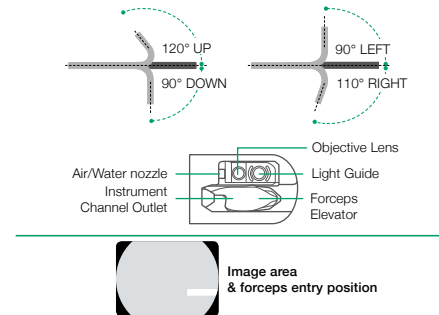

SUPER CCD TECHNOLOGY
The Super CCD and high performance optical system ensures high quality images.



FICE TECHNOLOGY
FICE can enhance slight colour differences such as vascular and mucosal patterns without tissue staining. The procedure digitally selects three wavelengths of light and displays reconstructed images.



HD TECHNOLOGY
Combine equipment displaying this logo to ensure that you view HDTV images on your monitor.

ACCESSORIES



SB-605


Suction Valve

AW-603

Air/Water Valve

AW-604G

Gas/Water Valve

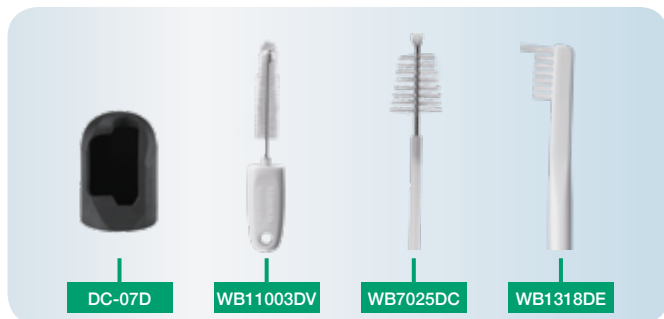

DC-07D

Distal End Cap

WB11003DV

Disposable Cleaning Brushes

WB7025DC

Disposable Cleaning Brushes

WB1318DE

Disposable Cleaning Brushes

**NEVER
STOP**



DUOMED

Duomed Swiss AG

Grenzstrasse 5a, 6214 Schenkon

Tel. +41 (0)41 510 07 00

www.duomed.com

FUJIFILM

FUJIFILM Europe GmbH

Heesenstr. 31, 40549 Düsseldorf, Germany

Tel.: +49 211-50 89 0, Fax: +49 211-50 89 8700

www.fujifilm.eu, endoscopy_eu@fujifilm.com