

Tortuga[®]


STONE EXTRACTION BALLOON



CLEAR STONES WITH CONFIDENCE

TORTUGA® features user-centric designs to ensure efficient stone extraction during ERCP. The versatile balloon accommodates both long and short guidewire techniques within a single balloon. The wider injection lumen enables effortless contrast injection.

The ergonomic handle is designed to offer an optimal grip and enhance user comfort. TORTUGA® is specifically engineered to maintain its robustness throughout the entire stone extraction procedure.

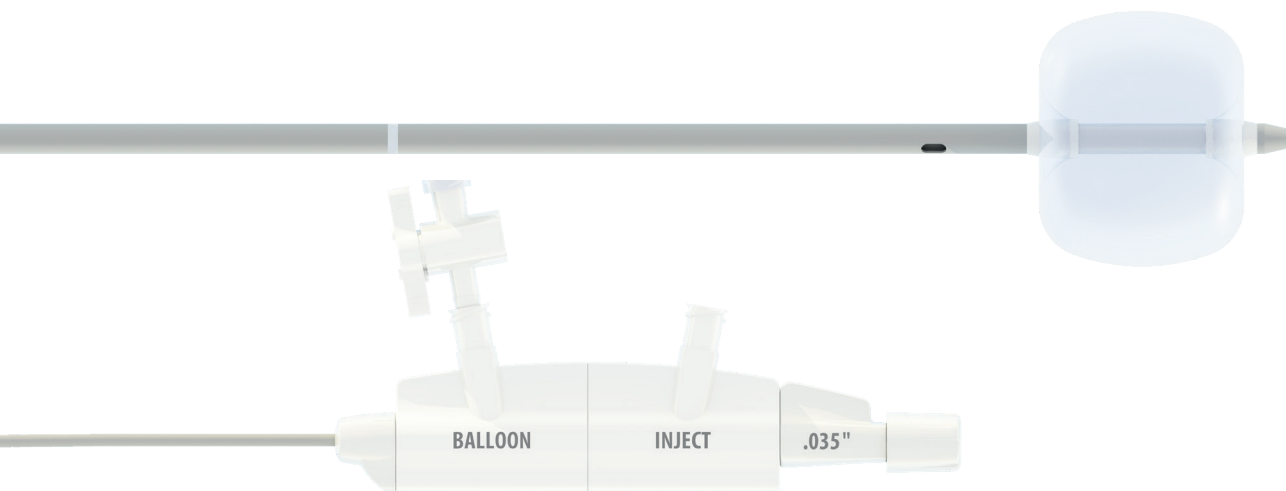

EXCELLENT BALLOON ROBUSTNESS

TORTUGA® is specifically engineered for high balloon robustness and has undergone extensive testing to ensure balloon burst resistance and uncompromising durability.

PRECISE 3-STAGE BALLOON INFLATION

Each balloon can be seamlessly adjusted to three different sizes with a single syringe:


- 8.5-12-15 mm
- 15-18-20 mm

ERGONOMIC ONE-PIECE HANDLE WITH CLEAR MARKINGS

The handle of TORTUGA® is designed to enhance the user's comfort during the procedure and its ergonomic shape fits naturally in the hand providing optimal grip. The clear distinction of ports and indicating markers supports an accurate and efficient procedure.

CLEAR FLUOROSCOPIC VISIBILITY



Two distinct X-ray markers at the distal and proximal ends of the balloon clearly indicate the balloon position on fluoroscopic view.

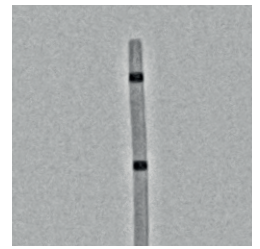




EFFORTLESS CONTRAST INJECTION

The wide contrast lumen enables smooth and easy contrast flow with less force.

(Left: TORTUGA®, Right: Conventional balloon)

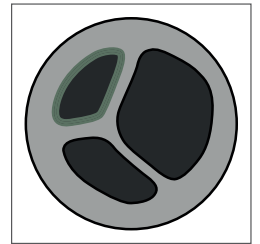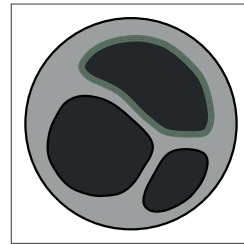


ONE BALLOON FOR BOTH SHORT AND LONG GUIDEWIRE

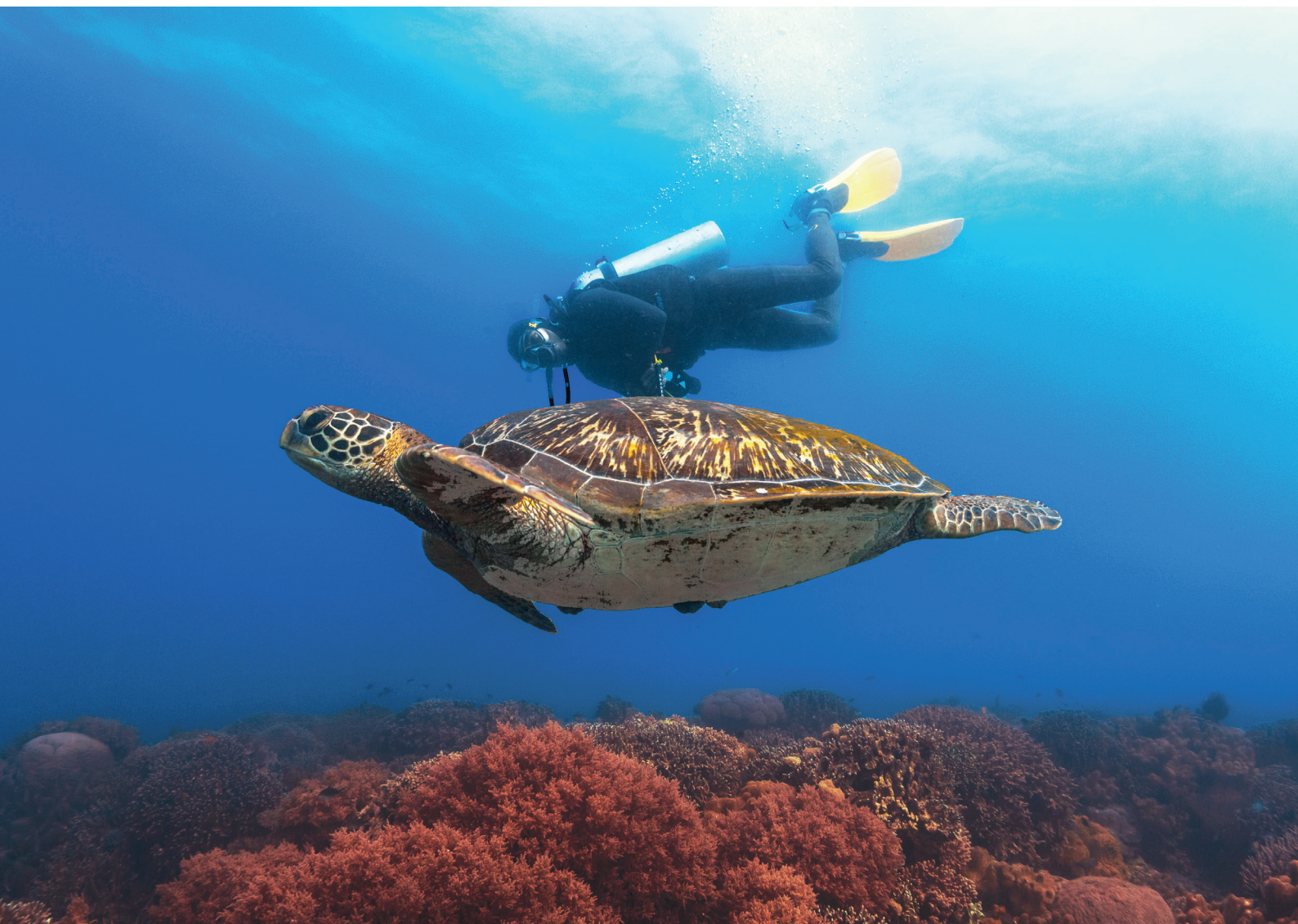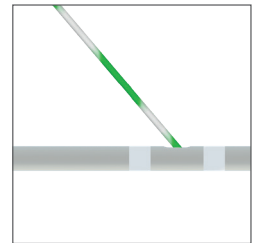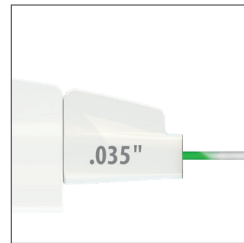


Product code	Material code	Injection	Balloon diameter (mm)	Length (mm)	Ø Working channel (mm)	Compatible guidewire	Unit
BAL1-D3-15-35-200	70100166445	Above	8.5-12-15	2200	3.8	0.035"	5
BAL1-D3-20-35-200	70100166446	Above	15-18-20	2200	3.8	0.035"	5
BAL1-P3-15-35-200	70100166449	Below	8.5-12-15	2200	3.8	0.035"	5
BAL1-P3-20-35-200	70100166447	Below	15-18-20	2200	3.8	0.035"	5

INJECTION ABOVE THE BALLOON

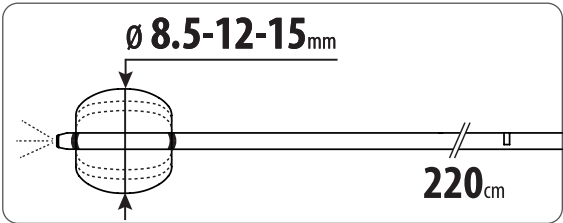


BAL1-D3-15-35-200

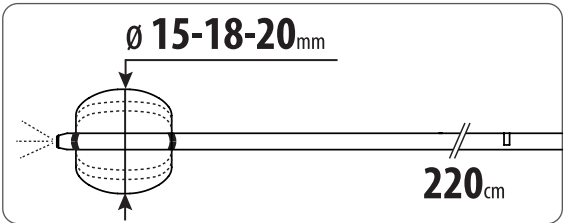


BAL1-D3-20-35-200

INJECTION BELOW THE BALLOON

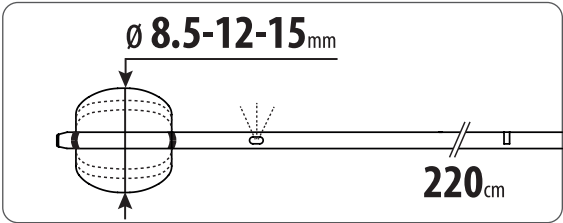

BAL1-P3-15-35-200

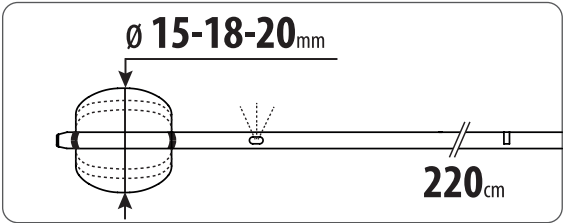

BAL1-P3-20-35-200



DUOMED
by Paalex

Duomed Swiss AG
Grenzstrasse 5a, 6214 Schenkon
Tel. +41 (0)41 510 07 00
www.duomed.com

FUJIFILM
Value from Innovation

Manufactured and distributed by
FUJIFILM medwork GmbH
Medworkring 1, 91315 Höchstadt, Germany
www.medwork.com

Sales Contact:
FUJIFILM Healthcare Europe GmbH
Balcke-Dürr-Allee 6, 40882 Ratingen, Germany
endoscopy_eu@fujifilm.com

© 2025 FUJIFILM Healthcare Europe

Specifications are subject to change without notice.
All brand names or trademarks are the property of their respective owners.
In some countries, regulatory approval may be required to import medical devices.
For proper use of the products, be sure to read the instruction for use prior to use.
For the availability of these products, please contact your local sales representatives.