


 Claritron®

# Single-use Ureteroscope

***SUV-4C-B 7.5Fr with Suction***

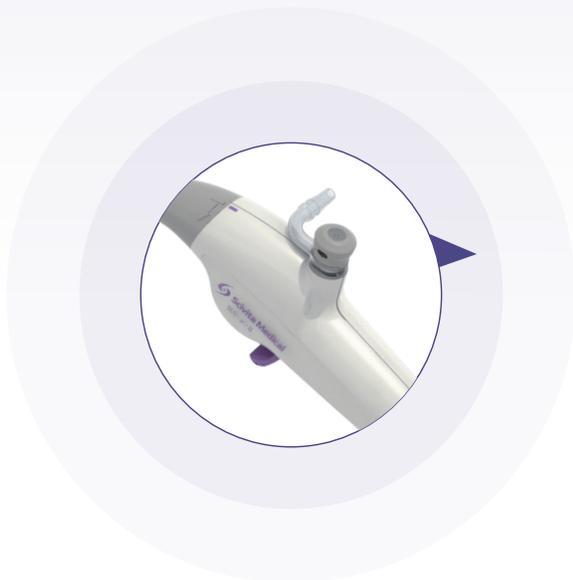

# Suction Liquid

## ● Superior Suction

It effectively clears liquid from renal pelvis, ureter and other areas, significantly enhancing the efficiency and safety of surgeries.

## ● Precise Control

The device allows for accurate control of the irrigation suction flow, thereby improving surgical precision and reducing risks.

## ● Enhanced Visual Field&Reduced Intra-renal Pressure

Through irrigation and suction, debris and blood clots are effectively removed, optimizing the visual field and minimizing intrarenal pressure.

Product Code	SUV-4C-B
Outer Diameter of Insertion Section	7.5Fr/Φ2.5mm
Outer Diameter of Distal End	7.2Fr/Φ2.4mm
Working Length	680mm
Total Length	990mm
Working Channel (ID)	3.6Fr/Φ1.2mm
Field of View	120°
Depth of Field	0.5~120mm
Deflection	285° Up/Down



SCIVITA MEDICAL TECHNOLOGY CO.,LTD.

☎ +86-512-8187 7788  
 🌐 www.scivitamedical.com  
 ✉ info@scivitamedical.com

No.2, Qingqiu Street,  
 Suzhou Industrial Park,  
 Suzhou, Jiangsu Province,  
 P.R.China, 215127

