


## Astopad® SOF Warming Cushion

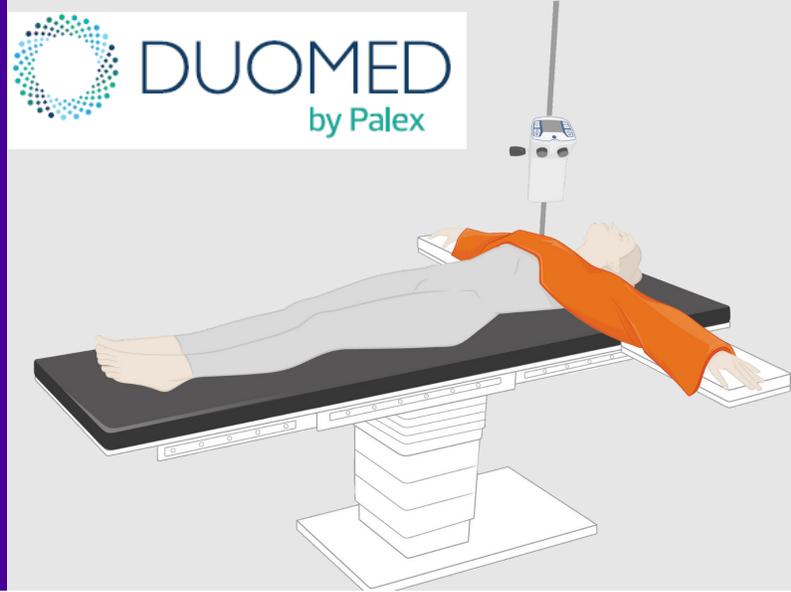


### OR TABLE PAD WITH INTEGRATED HEATING ELEMENT

The viscoelastic foam is combined with a dynamic foam base for ultimate decubitus prophylaxis on the operating table.

The integrated heating element on top of the OR table pad provides perfect warmth for the patients to prevent hypothermia.

Ideally, new operating tables can be directly equipped with Astopad SOF applied parts for the backrest pad or the entire surface of the operating table.

Astopad SOF pads are produced individually according to customer needs to fit perfectly within the size specifications of operating tables.

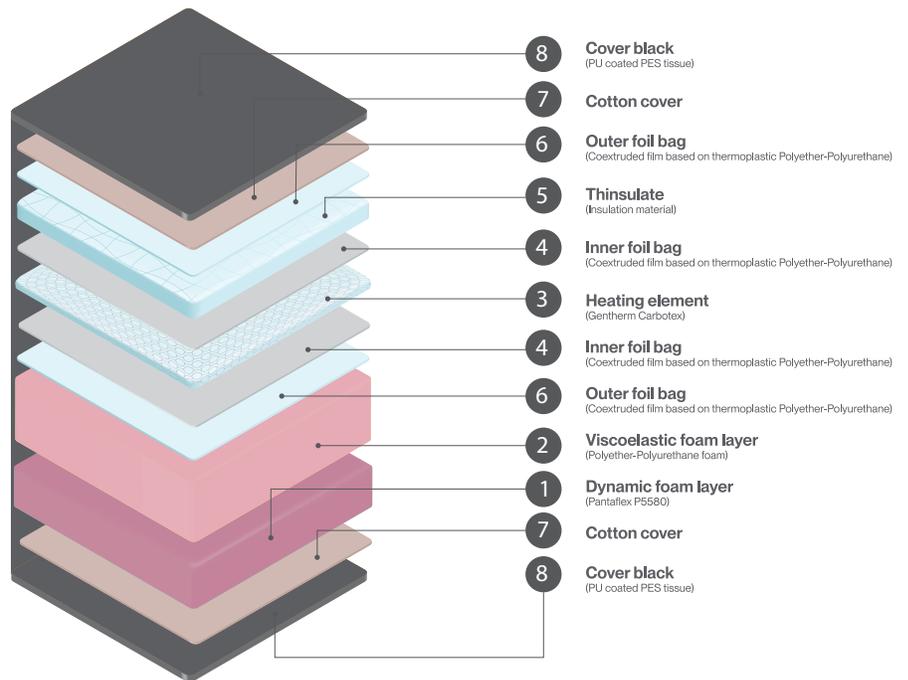

Patient Warming with Sustainability. Astopad applied parts are reusable, significantly reducing hospital waste.

### Features and Benefits



Safe and reliable design with 8 temperature sensors.



The Astopad system does not interfere with HF Surgical instruments.



Custom made to fit most OR tables.



Best possible decubitus prophylaxis with foam height up to 10 cm.



Different sizes available to fit all market relevant OR tables.



Silent operation for minimizing distractions in the OR.



Available as a backrest-cushion or mattress for the entire table.



Modular use together with COV possible.

## PRODUCT SPECIFICATIONS

**Controller Dimensions:**  
300mm H x 155mm W x 130mm D

**Battery Dimensions:**  
22mm H x 77mm W x 150mm D

**Weight (without battery):**  
2.0kg

**Weight (with battery):**  
2.5kg

### HEATING SYSTEM

**Heater Power:**  
160 Watts

### ELECTRICAL SYSTEM

**Electrical Connection:**  
100-240 VAC; 50-60 Hz 24 VDC

**Rated Current:**  
110V = 1.6 Amps / 240V = .8 Amps

**Heater Power:**  
160 Watts

**Time to Heat:**  
23°C to 37°C in approximately  
10 minutes

**Power Cord:**  
Approximately 250cm PVC Cable

# Astopad® System Ordering Information

CAT #	Product Description
<b>RESISTIVE WARMING DEVICE</b>	
DUO310	Astopad Controller (100V, 240V, 50-60Hz)
<b>RESISTIVE WARMING CUSHIONS</b>	
SOF2	Astopad Warming Cushion (mattress length) (2040-2300mm x 500-600mm x 40-100mm)**
SOF7	Astopad Warming Cushion (680-900mm x 440-600mm x 40-100mm)*
SOF4	Astopad Warming Cushion (1030-1300mm x 480-600mm x 40-100mm)*
SOF5	Astopad Warming Cushion (1480-1700mm x 480-660mm x 40-100mm)**
<b>ACCESSORIES</b>	
SOF402	Reusable cover for SOF2
SOF407	Reusable cover for SOF7
SOF404	Reusable cover for SOF4
SOF405	Reusable cover for SOF5

\* Heated, pressure relieving backrest for operating tables, antistatic. The dimensions are adapted individually, so that it is suitable for all operating tables.

\*\* Heated, pressure relieving operating table mattress, antistatic. The dimensions are individually adapted to the operating table.



Distributed by



T +45 44 47 47 48 • info.dk@duomed.com  
[www.duomed.com](http://www.duomed.com) • [in](#) Duomed Scandinavia A/S