



CLUTCHCUTTER



THE 3-IN-1 ESD TOOL FOR VARIOUS TECHNIQUES

The ClutchCutter is a well proven endoscopic mucosal dissection tool that can enable the user to grasp, and electrically dissect and coagulate the mucosa of the upper and lower GI tract. In procedures that need to be performed perpendicular to the muscular layer, the ClutchCutter makes it possible to grasp, lift up and incrementally dissect the tissue.

ClutchCutter



**TOOTHED JAWS TO GRIP
THE MUCOSA EFFICIENTLY**


**ROTATABLE DISTAL JAWS FOR
A PRECISE LESION APPROACH**

**INSULATED OUTER EDGE TO REDUCE
RISK OF DAMAGING LOCAL TISSUE**


CLUTCH CUTTER SINGLE USE ESD KNIFE



Model	Material Code	Characteristics	Identifier	Knife length (mm)	Working length (mm)	Ø Working channel (mm)
DP2618DT - 35 -	16444210	ESD Scissor	- 35 -	3.5	1,800	2.8 or more
DP2618DT - 50 -	16444222	ESD Scissor	- 50 -	5.0	1,800	2.8 or more

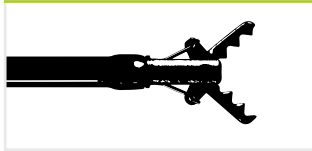

DP2618DT -35-: 3.5 mm jaw length

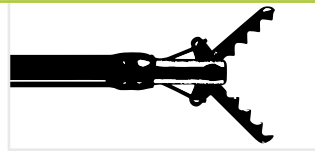

DP2618DT -50-: 5.0 mm jaw length

3-in-1 ability

All stages of ESD can be performed with this 3-in-1 tool, from marking to making incisions and coagulating.



1. Incision



2. Dissection



3. Coagulation